


2'-O-Aminolinker-5-methyl U CEP
Product No. BA 0191

Product Information

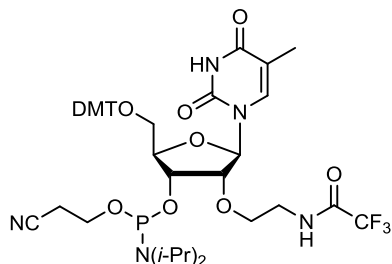

$C_{44}H_{53}F_3N_5O_{10}P$
Mol. Wt.: 899.89

Coupling, cleavage, and deprotection: 2'-O-Aminolinker-5-methyl U CEP couples with greater than 95% efficiency (typically >97%) using the standard RNA protocols recommended for popular synthesizers. Cleavage and deprotection may be accomplished using standard techniques e.g., concentrated ammonium hydroxide at 55 °C for 18 h.