


2'3'-Di-O-acetyl-U-5'-CEP

Product No. BA 0358

Product Information

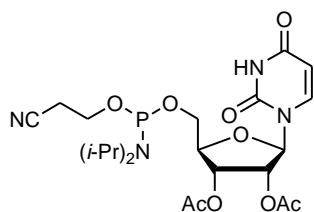

C₂₂H₃₃N₄O₉P
Mol. Wt.: 528.49

Use: Dissolve the phosphoramidite in acetonitrile at concentrations recommended by the synthesizer manufacturer. Coupling should be carried out using standard instrument RNA protocols using a 12 minute coupling time. Cleavage from the solid support can be carried out under standard conditions, and standard deprotection conditions may be employed. Coupling efficiencies of > 95% can be obtained with BA 0358.